

One device. Diverse patients. Customised solutions.


APA™ is an all-in-one traditional and Video Laryngoscope for routine, paediatric, rescue and difficult airways. Its multiple interchangeable components permit you to switch between MAC, MIL, Oxygenation and Difficult Airway allowing you to build the airway device you require.




Wide display



Anti-glare

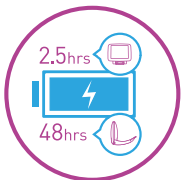

Dual battery system

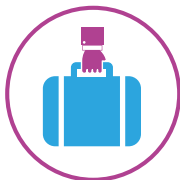


Portable



Diverse groups



Continuous oxygen delivery

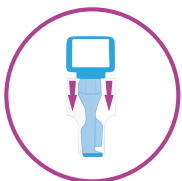

Increased protection

MIL – Paediatric

Introducing the new Miller style interchangeable Camera Module and blades to help your paediatric intubations. Paediatric blades are available in sizes MIL 1 and MIL 2.



MAC – Adult

APA™ can be used as both a traditional laryngoscope and as a Video Laryngoscope for improved visualisation to assist you in securing the airway. For your routine airways a MAC 3 and MAC 4 blade are available.



Difficult Airway

A unique prism directs the line-of-sight around the anatomy. The blades are designed to assist in lifting the epiglottis for Cormack-Lehane grade III and IV cases. Our original DAB (Difficult Airway Blade) offers a built-in tracheal tube channel to assist the ETT into the laryngeal opening and we also offer the U-DAB (Unchannelled Difficult Airway Blade).



Oxygenation

MAC 3 and MAC 4 oxygenation blades assist routine endotracheal intubation by delivering continuous oxygenation during laryngoscopy, when additional apnoeic oxygenation is clinically advantageous.



Product code	Description	Qty/Box
Video Laryngoscope		
701000	APA™ Video Viewer and Charger	1
701006	APA™ Handle	1
700930	APA™ MAC Camera Module	1
702910	APA™ MIL Camera Module (S)	1
702930	APA™ MIL Camera Module (L)	1
Disposable Blades		
702110	APA™ MIL 1 Blade	10
702120	APA™ MIL 2 Blade	10
700130	APA™ MAC 3 Blade	10
700140	APA™ MAC 4 Blade	10
700170	APA™ DAB	10
700180	APA™ U-DAB	10
700330	APA™ O ₂ MAC 3 Blade	5
700340	APA™ O ₂ MAC 4 Blade	5
Accessories		
701002	APA IP Shield™	50
701004	APA™ Video Cable (3.0m)	1
701008	APA™ Carry Case	1

All blades are clean, non-sterile and anti-fog coated.

